

# Data Sheet

# High Capacity Manual Force Test Stands TSF & TSFH

32-1055 REV 3

The TSF test stand is ideal for spring testing, destructive testing, product integrity testing, and other applications requiring up to 1,000 lbf (5,000 N) of push or pull force. The stand features a large side mounted hand wheel requiring minimal user effort and a rigid 3" x 3" column. A removable steel loading surface with threaded holes for fixture and grip mounting adds flexibility for many testing applications, and a solid metal base (TSF) contains a storage pocket and holes for bench mounting. Available in vertical and horizontal orientations.

## Specifications

Load Capacity:	1000 lbf [5000 N]
Maximum Travel:	4.0 in [102 mm]
Loading method:	Side-mounted hand wheel
Travel rate:	0.013 in [0.34 mm] per wheel rev.
Weight:	TSF: 33 lb [14.5 kg] TSFH: 27 lb [12.3 kg]
Included accessories:	Tool kit, three hooks (small, medium, and large), 2" and 3" compression plates, two couplers (#10-32 and 5/16-18), force gauge mounting hardware
Warranty:	3 years (see individual statement for further details)

## Optional Equipment


Digital travel display TSF001	6 in [150 mm] travel, 0.0005 in [0.01 mm] resolution, 0.001 in [0.03 mm] accuracy, SPC output for data collection.
Horizontal / wall mounting kit TSF002	Converts the TSF to a TSFH
Indicator mounting kit AC1004	If a indicator with remote sensor is used with the TSF/TSFH, the indicator can be securely mounted to the side of the column with this adapter.

» Digital travel display



## Dimensions


in [mm]

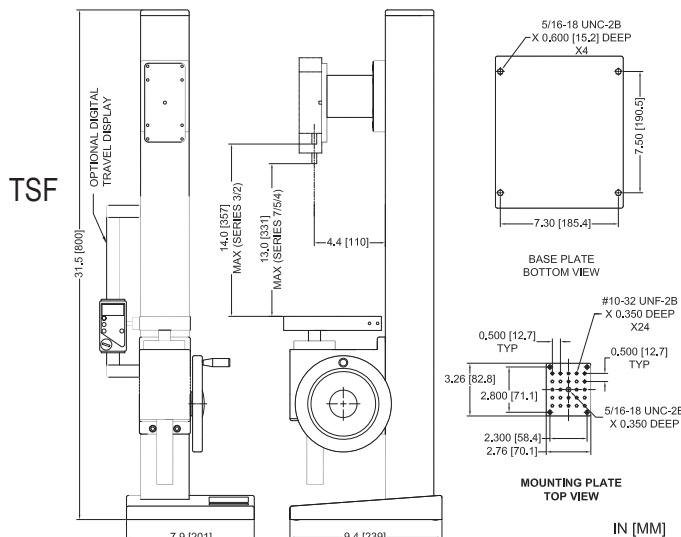

TSF




TSFH

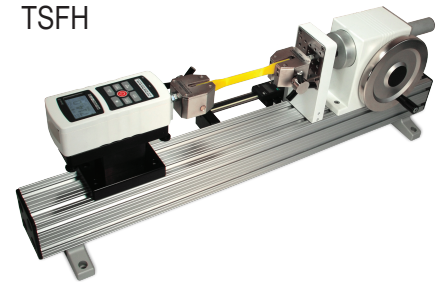

The TSF is shown testing a compression spring, with a Series 5 digital force gauge and digital travel display

The TSFH is shown in a tensile test, with Series 5 digital force gauge, film and paper grips, and digital travel display

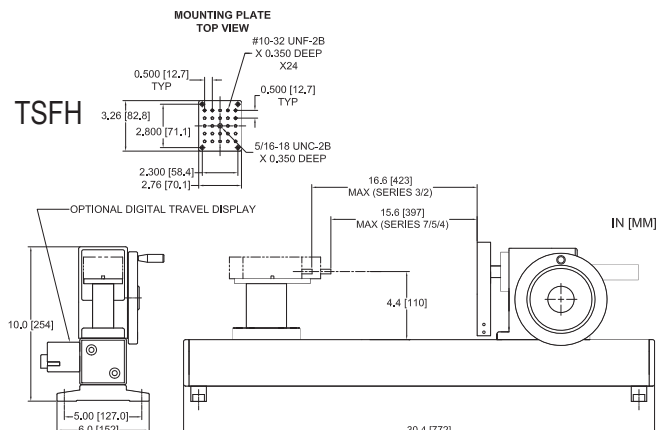

## Ordering Information

Model No.	Description
TSF	Test stand, hand wheel-operated, 1000 lbf, vertical
TSFH	Test stand, hand wheel-operated, 1000 lbf, vertical
TSF001	Digital travel display for TSF/TSFH
TSF002	Horizontal / wall mounting kit for TSF/TSFH
AC1004	Mounting kit, Indicator to test stand



Mark-10 Corporation ■ www.mark-10.com ■ info@mark-10.com  
Toll-free: 888-MARK-TEN ■ Tel: 631-842-9200 ■ Fax: 631-842-9201



Specifications subject to change without prior notice