



EYE-LEVEL DIGITAL CLINICAL SCALES

- 600 lb x 0.2 lb / 300 kg x 0.1 kg capacity
- Sonar touchless or inline mechanical height rod
- Extra-wide and flat platform (17 in W x 17 in D x 2.75 in H)
- Wi-Fi/Bluetooth models available
- Units locking into kilograms or pounds



WEIGHT + HEIGHT INSTANT READOUTS WITHOUT PUSHING A SINGLE BUTTON

Sonar touchless height rods available on some apex® models provide immediate digital readouts for optimal speed, accuracy, and hygiene

Model APEX-SH

Model APEX



OF PATIENT
MEASUREMENT
EFFICIENCY
+ STYLISH
DESIGN
IS NOW AVAILABLE!

STATE-OF-THE-ART SONAR HEIGHT RODS

DETECTO's unique touchless sonar height rod instantly locks on the patient's height and measures nearly any patient size. The apex's sonar height rod measures from 24 inches (2 ft) up to 87.5 inches (7'3") by 0.2 inch increments and metric readouts of 60 cm up to 222 cm by 0.5 cm increments (model APEX-SH). Simply step on the scale and your height will be captured and displayed on the LCD within a second.




DETECTO's apex®
scales measure weight
and height efficiently
on patients of all
ages and sizes.



The apex's inline mechanical height rod models measure from 24 inches (2 ft) up to 87.5 inches (7'3") by 0.125 inch increments and metric readings of 60 cm up to 222 cm by 0.1 cm increments (model APEX). To calculate BMI, simply use the scroll arrows to enter the height measurement and then press the BMI key.

- Eye-level weight display convenience
- 10 different standard apex® models available for you to choose from
- Sturdy column construction



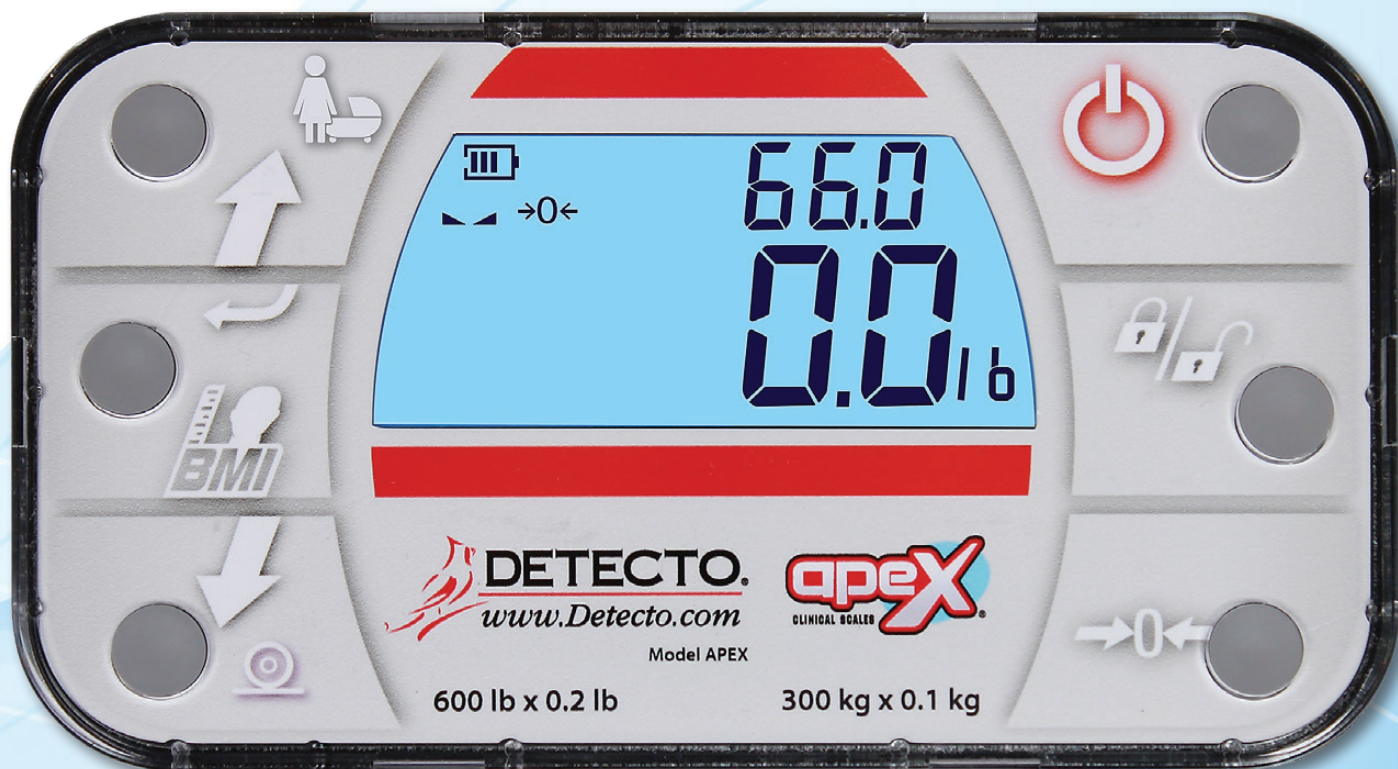

HIGH-TECH WEIGHT INDICATOR


The apex® features a high-tech, smart-phone-style indicator for easy utilization with 6 simple buttons and a bold LCD display for large readouts.

- 0.75-inch-high, clinical-blue LCD weight readouts
- Body Mass Index calculation
- Weight and height displayed onscreen simultaneously
- Automatically locks on weight temporarily for a few seconds in order to easily view and record measurement
- Battery power level indication
- Mother/Baby function factors the weight of infants and toddlers held by an adult



Wi-Fi/Bluetooth models available

Dual integral wheels provide for easy transport on all apex® models.

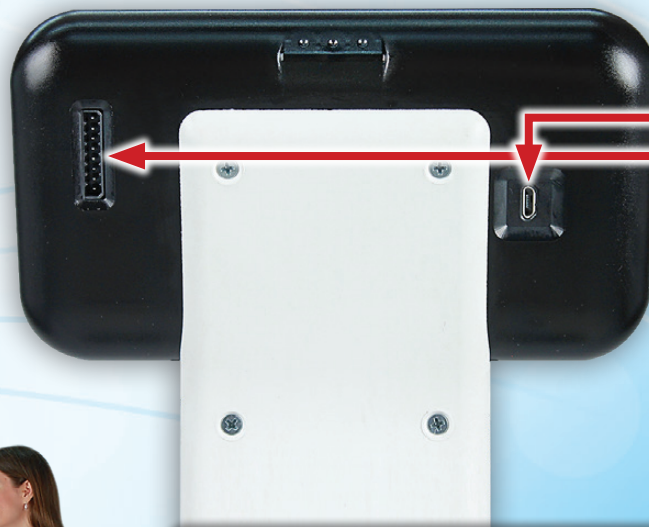

VERSATILE CONNECTIVITY:

- 1 micro USB-B port
- 1 RS232 serial port
- HL7 IEEE 11073 compliant (standard)
- Wi-Fi/Bluetooth models available for wireless EMR/EHR
- Welch Allyn® device connectivity available: CVSM, CSM, and LXI (see models chart on back page)

POWER



apex® models are available with the AC adapter included standard (see models chart on back page).

An integral bubble level helps locate the most level surface area on your floor.

The lift-off apex® platform cover may be removed for wipe-down cleaning or to replace the 12 AA batteries.

The scale's structural base integrity is manufactured rock-solid for optimal patient safety and stability.



apex® EYE-LEVEL DIGITAL SCALES WITH MECHANICAL HEIGHT RODS						
Model		APEX	APEX-AC	APEX-C-AC	APEX-UWA-AC	APEX-LXI-AC
Capacity		600 lb x 0.2 lb / 300 kg x 0.1 kg				
Height Rod		Mechanical inline height rod				
Height Rod Measures		24" (2') - 87.5" (7'3") x 0.125 in / 60 cm - 222 cm x 0.1 cm				
DIMENSIONS	Overall Dimensions	17 in W x 19.5 in D x 53.3 in H / 43 cm W x 49.5 cm D x 135.3 cm H				
	Platform Size	17 in W x 17 in D x 2.75 in H / 43.2 cm W x 43.2 cm D x 7 cm H				
	LCD	.75 in / 1.9 cm H LCD digits				
POWER	Batteries	12 AA batteries (not included)				
	AC Adapter (included)					
Ports		1 RS232 serial port and 1 micro USB-B port				
Connectivity				Wi-Fi / Bluetooth	Welch Allyn CVSM/CSM	Welch Allyn LXI
Net Weight		31 lb / 14 kg				
Shipping Weight		40 lb / 18 kg				
Shipping Dimensions		58 in x 20 in x 6 in / 147 cm x 51 cm x 15 cm				
UPC Code		809161197808	809161197907	809161322606	809161198102	809161198201

apex® EYE-LEVEL DIGITAL SCALES WITH SONAR HEIGHT RODS						
Model		APEX-SH	APEX-SH-AC	APEX-SH-C-AC	APEX-SH-UWA-AC	APEX-SH-LXI-AC
Capacity		600 lb x 0.2 lb / 300 kg x 0.1 kg				
Height Rod		Touchless sonar height rod				
Height Rod Measures		24" (2') - 8.5" (7'3") x 0.2 in / 60 cm - 222 cm x 0.5 cm				
DIMENSIONS	Overall Dimensions	17 in W x 19.5 in D x 89.5 in H / 43 cm W x 49.5 cm D x 227.3 cm H				
	Platform Size	17 in W x 17 in D x 2.75 in H / 43.2 cm W x 43.2 cm D x 7 cm H				
	LCD	.75 in / 1.9 cm H LCD digits				
POWER	Batteries	12 AA batteries (not included)				
	AC Adapter (included)					
Ports		1 RS232 serial port and 1 micro USB-B port				
Connectivity				Wi-Fi / Bluetooth	Welch Allyn CVSM/CSM	Welch Allyn LXI
Net Weight		32 lb / 14 kg				
Shipping Weight		41 lb / 19 kg				
Shipping Dimensions		58 in x 20 in x 6 in / 147 cm x 51 cm x 15 cm				
UPC Code		809161197808	809161197907	809161322804	809161198102	809161198201

DETECTO reserves the right to improve, enhance, or modify features and specifications without prior notice.



102 E. Daugherty, Webb City, MO 64870 USA
Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-2153

SOLD BY: